**[ comfort zone ]**

**MICROSONIC  
FACIAL DEVICE**

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### 1 BATTERY CHARGER

### 1 USER MANUAL

PLEASE NOTE: read the instructions in the manual carefully before using the device as it contains important information regarding safety, usage, maintenance and warnings.

Store the manual carefully in case you need to check back in the future.

This device should only be used for the purpose it was designed for, i.e. as a microcurrent stimulator for the face and an ultrasound massager for aesthetic purposes and personal use. It should only be used according to the methods set out in this manual.

Any other forms of use are inappropriate. The manufacturer cannot be held responsible for any damages deriving from improper, incorrect or unreasonable use of the device.

## MICROSONIC FACIAL DEVICE


An exclusive skincare device with double ULTRASOUND and MICROCURRENT technology, used by professionals, with three main effects:

- 1. deep cleansing of the skin**, to remove the top layers of the skin and prepare the skin for subsequent treatments. (SONIC mode)
- 2. regenerative anti-aging stimulation** of the deeper layers of the skin using microcurrents, which serve as a powerful stimulator for cellular metabolism and cutaneous circulation. (MICRO mode)
- 3. optimization of product absorption and regenerative massage** which significantly improves skin hydration and facilitates the absorption of the active ingredients of any cosmetics applied. (SONIC + MICRO mode)

The device delivers a complete, highly effective facial treatment in the privacy and comfort of your own home. It offers a range of modes based on technology, intensity and the position of the blade.

MADE IN ITALY

## HOW TO USE THE DEVICE

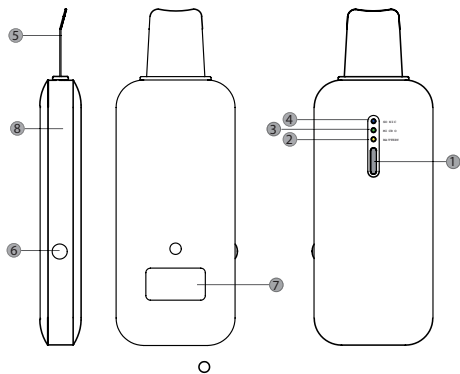1. Operating button
2. Charging LED
3. Microcurrent LED
4. Ultrasound LED
5. Application blade
6. Charging port
7. Data label
8. Metal border  
(microcurrent electrode)

### TURNING ON THE DEVICE


In order to turn the device on, hold down the button **(1)** for at least 3 seconds. A sound will notify you that the device will come on. The device will go into SONIC mode (blue LED on) **(4)**.

### SELECT DESIRED FUNCTION


By quickly pressing button **1**, you can bring up the menu and select your desired function from the following sequence:

**SONIC:**  
blue LED on  
**(4)**

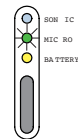**MICRO:**  
green LED on  
**(3)**

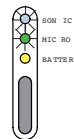
**SONIC+MICRO**  
blue and green  
LED lights on  
**(3+4)**

Once you have decided the mode you want to use, this will automatically be selected if you do not press the button for a further 5 seconds.

After this period, the button can be used to adjust the intensity or change the mode.

### POSITION OF THE BLADE ON THE FACE

Based on the program being used, the device will be held in position 1 (blade against skin) or position 2 (blade resting on skin).
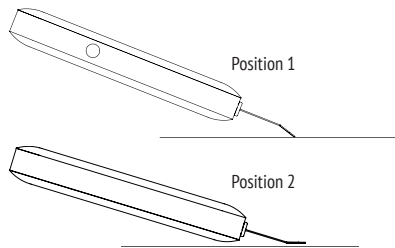
## TECHNOLOGIES

### SONIC TECHNOLOGY

#### image of blade in position 1

Through intense ultrasonic mechanical vibrations, this technology enables the device to have a deep cleansing and scrubbing effect. When the SONIC technology is activated, the BLUE LED light will come on (constant for normal and hard skin, flashing for sensitive skin).

For this mode, the blade is held in position 1. You must moisten the parts of your face that you are going to use the device on with a damp cotton pad or a water-based solution. Next, pass the blade over the desired areas of the face, exerting some light pressure. Make repetitive movements over the desired areas. You will notice that some of the water evaporates: this is a sign that the device is working correctly.

This facial treatment should last for 4 minutes.

### MICRO TECHNOLOGY

#### image of blade in position 2

When this technology is activated, the device begins emitting delicate yet effective microcurrents with alternative positive and negative impulses. These act on cellular metabolism and cutaneous circulation.

When the MICRO technology is activated, the green LED light will come on: fast flash for hard skin, slow flash for sensitive skin, medium flash for normal skin.

For this mode, the blade is held in position 2. You must moisten the parts of your face that you are going to use

the device on with a damp cotton pad or a water-based solution. Next, pass the blade over the desired areas of the face, exerting some light pressure. Make repetitive movements over the desired areas. During the treatment, you must make sure that your skin is well-moistened in order to ensure that there is good electrical contact between the blade and the skin. You should wet your face several times during the process. The metal border around the edge of the device is the electrode for the microcurrent. This treatment should last for 4 minutes.

### SONIC + MICRO TECHNOLOGY

#### image of blade in position 2

Combining the benefits of both technologies (microcurrents and ultrasound) is the perfect way to optimize absorption of cosmetic products that are applied, to regenerate the skin and to improve the tone of the skin or to help apply a [ comfort zone ] product.

When the SONIC + MICRO technologies are activated, the GREEN and BLUE LED lights are both on.

During usage in MICRO and SONIC + MICRO modes, hold the blade in position 2, apply your selected product and pass the blade over the area until absorbed. Ensure that you grip the device well, making sure that you maintain good contact between your hand and the metal border, in order to ensure that the microcurrent is being passed through.

This facial treatment should last for 5 minutes.

# TREATMENTS

## SENSITIVE SKIN

### 1. Cleansing and scrubbing

Hold down the button for 3 seconds to turn on. The blue LED will come on. Wait for the 2 beeps, then press again (intermittent blue LED).

Wet the skin with water, put the blade in position 1 and start the treatment. If necessary, wet the skin again.

### 2. Regenerative massage

Hold down the button for 3 seconds to turn on. Press once: the green LED will come on (slow flash).

Wet the skin with water, put the blade in position 2 and start the treatment. If necessary, wet the skin again.

### 3. Product absorption and regenerative massage

Hold down the button for 3 seconds to turn on, then press it twice: the green and blue LED lights will come on.

Apply your chosen product, put the blade in position 2 and start the treatment.

## NORMAL SKIN

### 1. Cleansing and scrubbing

Hold down the button for 3 seconds to turn on. The blue LED will come on and stay on.

Wet the skin with water, put the blade in position 1 and start the treatment. If necessary, wet the skin again.

### 2. Regenerative massage

Hold down the button for 3 seconds to turn on. Press once: the green LED will come on. Wait for the 2 beeps, then press again. The green LED will come on (medium flash). Wet the skin with water, put the blade in position

2 and start the treatment. If necessary, wet the skin again.

### 3. Product absorption and regenerative massage

Hold down the button for 3 seconds to turn on, then press it twice: the green and blue LED lights will come on.

Apply your chosen product. Put the blade in position 2 and start the treatment.

## HARD SKIN

### 1. Cleansing and scrubbing

Hold down the button for 3 seconds to turn on. The blue LED will come on. Wet the skin with water, put the blade in position 1 and start the treatment. If necessary, wet the skin again.

### 2. Regenerative massage


Hold down the button for 3 seconds to turn on. Press once: the green LED will come on. Wait for the 2 beeps, then press twice. Wet the skin with water, put the blade in position 2 and start the treatment. If necessary, wet the skin again.

### 3. Product absorption and regenerative massage

Hold down the button for 3 seconds to turn on, then press it twice: the green and blue LED lights will come on.


Apply your chosen product, put the blade in position 2 and start the treatment.

## SPECIFIC TREATMENTS

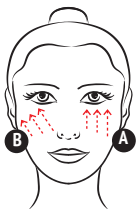
### CHEEK

**A** vertical movements from bottom to top over the cheekbone  
**B** horizontal movements, from the edge of the lip to the ear, as if gently "lifting" the skin.

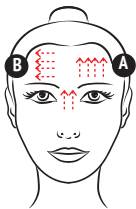
### FOREHEAD AREA

The blade can be moved slowly in either of the modes over the forehead and between the eyebrows:  
**A** vertical movements from bottom to top  
**B** horizontally from the temples to the outer part of the forehead.

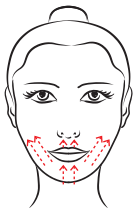
### JAW, CHIN AND UPPER LIP AREA

You must gently move the blade from bottom to top, from the chin to the ear, following the arch of the jawbone, from the tip of the chin to the bottom lip and from the upper lip to the start of the nose.
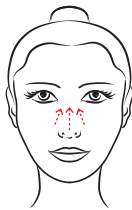
### NOSE

The deep cleansing effect made possible by the device is particularly important for the entire nose area, which often attracts blackheads. These can be removed if done carefully, while effective exfoliation can help prevent grease and oil. The ultrasound micromassage and the application of an appropriate cosmetic product complete the treatment. The microcurrent regeneration treatment is not usually necessary in this area.
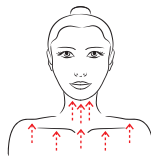
### NECK AND CHEST AREA

This is a particularly delicate area. Rather than cleansing and peeling (which are still possible), micro-massage and microcurrent tissue regeneration are the most effective treatments here. It's important to make gentle movements from bottom to top.

# TECHNICAL INSTRUCTIONS

## CLEANING THE DEVICE

Use a soft cloth to clean the device, adding a spray of neutral detergent to the cloth if required. Do not use alcohol, diluents, solvents or other chemical products.

## CHARGING THE BATTERY

When the battery life dips below a critical level during usage, the yellow LED (2) will start flashing to indicate that the battery needs recharging. If you continue to use the device, it will eventually turn off automatically.

When you insert the charger into the port, (6) the yellow LED (2) will come on and remain on until the battery is fully charged. It takes around 4 hours for the battery to charge fully.

During charging, the device cannot be used for safety reasons. Do not recharge the device while it is exposed to direct sunlight, heat sources or steamy environments.

Only used the included **VATSPA05** charger. Both the device and the charger will get slightly warmer during the charging process. This is normal and not a sign of malfunctioning.


Keep the charger away from water and only handle it with dry hands.

## REMOVING THE BATTERY

The device is fitted with a lithium-ion (Li-Po) polymer battery in the handle. The battery is rechargeable and – if used correctly – can last for over 300 charging/ discharging cycles.

The battery cannot be removed by the user. If the battery needs replacing, you must contact a Vupiesse-authorized center.

The procedure set out below is only to be used to remove the battery so that the device can be properly disposed off after its product life is over.
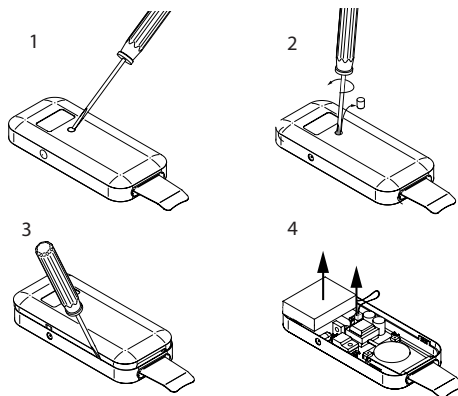
## DISPOSAL

The product is an electronic device and as such contains electronic components which may be harmful to the environment and to human health if disposed of in a way that doesn't comply with the law. Ensure that you dispose of the device itself, its component parts and any accessories in accordance with local waste disposal regulations. Make sure that you dispose of batteries in the special disposal containers once no longer of use. Directive 2012/19/EU obliges product owners to dispose of this product by taking it to their local Electrical and Electronic Equipment Waste center (WEEE) or by handing it in to the seller when they make a new purchase. The directive obliges product owners to ensure that they do not dispose of the device in the general waste. The directive sets out sanctions regarding improper disposal of this type of waste. The symbol below (wheelie bin with a cross through it) on the data label of the device indicates that electrical and electronic waste must be disposed of separately.

## CLASSIFICATION

- Device with internal power source; Class II when connected to charger
- Type BF applied section
- Grade of protection against penetration of liquids: IPX0

## PRECAUTIONS

Keep the skin moist using water, a specific water-based solution or a suitable [ comfort zone ] product during microcurrent treatments.

Do not immerse the device in water.

Keep out of reach of children.

## DO NOT USE IF YOU HAVE ANY OF THE FOLLOWING

- Pathologies of the face
- A pacemaker
- Metal prostheses or plates
- Neurological pathologies, including epilepsy and multiple sclerosis
- Tinnitus
- Debilitating psychophysical conditions
- Skin bruising or abrasions
- Open wounds or mucous issues
- Conditions where your skin is swollen, exuding pus/ discharge or sore
- Large moles that are protruding and irregular in shape
- Sunburn
- External use only

Consult your doctor before using the device if you have had recent cosmetic procedures involving chemical substances/chemical peel.

If you have any doubts, please consult your doctor.

## TECHNICAL INFORMATION

| Ultrasound parameter           | Ultrasound   |
|--------------------------------|--|
| Ultrasound carrier frequency   | 23-26Khz   |
| Intermittent modulation        | 0.8 Hz   |
| Low Intensity Freq-Duty Cycle  | 50Hz - 30%   |
| High Intensity Freq-Duty Cycle | 60Hz - 60%   |
|                                | <b>Microcurrent</b>  |
| Max current intensity          | 500 µA (1KOhm current)   |
| Microcurrent carrier frequency | 8 Hz   |
| Low Intensity Duty Cycle       | 17%  |
| Medium Intensity Duty Cycle    | 34%  |
| High Intensity Duty Cycle      | 50%  |
| <b>General parameters</b>      |  |
| Power source                   | 7.4V Internal lithium-ion polymer battery  |
| Weight                         | 180 g  |
| Size                           | 146 x 51 x 21 mm   |
| Usage conditions               | Temperature: 15 - 35 °C Humidity: 20 - 65% RH  |
| Storage conditions             | Temperature: 10 - 50 °C Humidity: 10 - 85 % RH   |
| Charging                       | VUPIESSE mod. GPE060A-120020-Z (Vupiesse reorder code: VATSPA05)<br>Input: AC 100-240V- 50-60Hz 0.2A Output: DC 12V - 200mA 2.4W |

## REGULATORY COMPLIANCE

The device has the CE marking. The design and production of the device comply with the requirements set out in the UNI EN ISO 9001 and UNI EN ISO 13485 quality standards. With regards to product checks, the **EN60335-1 / EN60335-2-32** electrical safety standards have been applied. Electromagnetic compatibility compliance has been verified through the application of the **EN62233 / EN55014-1 / EN55014-2** standards.

1. Symbol indicating device compliance with European Directives
2. 12V battery charger (direct current)
3. Device with double insulation
4. Electrical and Electronic Equipment Waste WEEE: devices with this label must be disposed of separately to general waste.

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Corporation