

WELCOME GUIDE


* For safety issue, the product remains low capacity during shipping. Please fully charge the product before initial use!

SmartDent
— OPC-ET1 —

Please read this Guide carefully before using the product

Package Contents

Toothbrush Handle x 1
Toothbrush Replacement Head x 2
Toothbrush Holder x 1
Type-C Cable x 1
Welcome Guide

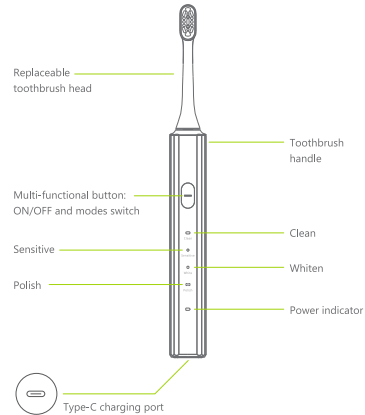

Functions

• Power On / Off

Press the **multi-functional button** to turn on/off the toothbrush.

💡 Built-in 2-minute timer will avoid brushing teeth for a long time. Toothbrush will also vibrate every 30 seconds to remind you to switch the cleaning area.

• Brushing Modes

The **LED light** indicates the brushing mode.

Clean: Effectively remove plaque. It is the standard mode.

Sensitive: Gently brushing for sensitive teeth and gums. This is the suggested mode for a new user.

White: Whiten the surface of the teeth.

Polish: It can not only whiten the teeth, but also massage gums and improve gum health.

• Mode Switch

Turn on the toothbrush, press the multi-functional button to switch the modes. Long press the multi-functional button more than 5 seconds will turn off the toothbrush.

💡 Memory feature of toothbrush means it will record your brushing habits and you won't have to reset the intensity level each time you brush.

• Charge Your Brush


Toothbrush Status Indicator		Status
● ○ ●	Slow red flashing	Low battery
● ○ ● ● ●	Rapid red flashing	Being charged
●	Steady red	Fully charged

💡 Please charge the tooth brush with 5V/1A charger(maximum 2A).

Steps for Using


Step 1

Fix the toothbrush head on the handle correctly.

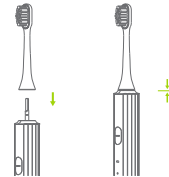


Step 2

Squeeze toothpaste on the bristle. Turn on the toothbrush after putting it into your mouth.

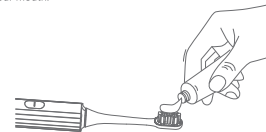


Step 3

Move slowly to brush different areas. Short vibration during brushing reminds you to change the cleaning area.

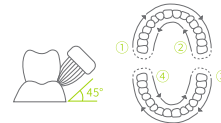

Step 4

Rinse the tooth brush head and bristle after each use.

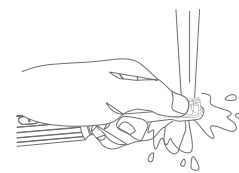

Specifications

Rated Input	5V=1A
Rated Voltage / Power	3.7V — /1.8W
Battery Capacity	600mAh
Charging Time	2 hours
Battery Life	20 days
Charging Type	Type-C
Brush Head Material	0.152mm Dupont Bristles
Dimensions	245.5mm*25.2mm*25.2mm

Attention

- Please make sure the port plug is tightly sealed when using.
- Toothbrush will stop working during charging.
- The brush holder must be attached at smooth, flat and clean wall or glass surface.

Important Safety Instructions

When using this product, basic precautions should always be followed, including the following:

1. Always place and store the toothbrush and charging cable at a safe distance from water.
2. Unplug the charging cable immediately if it falls into water.
3. Please make sure the port plug is tightly sealed.
4. Discontinue use if the product appears damaged in any way (brush head and handle).
5. Keep the toothbrush and cable away from heated surfaces.
6. Use this appliance only for its intended use as described in this document or as recommended by your dental professional.
7. This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
8. Children should be supervised to ensure that they do not play with the appliance.
9. Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
10. Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
11. Discontinue use of this product and contact a physician/ dentist if discomfort or pain is experienced.
12. Consult your physician if you have other medical concerns.