oraimo

WELCOME GUIDE

MagniSonic | OPC-ET500N |

*Please read the instructions carefully before using the product and keep the user manual in a safe place for future reference.

SCAN & DOWNLOAD

oraimo store

Easily Smart, Simply Best

Authentic Products

Fast Delivery


Best Price

Prompt Service

PACKAGE CONTENTS & PRODUCT OVERVIEW

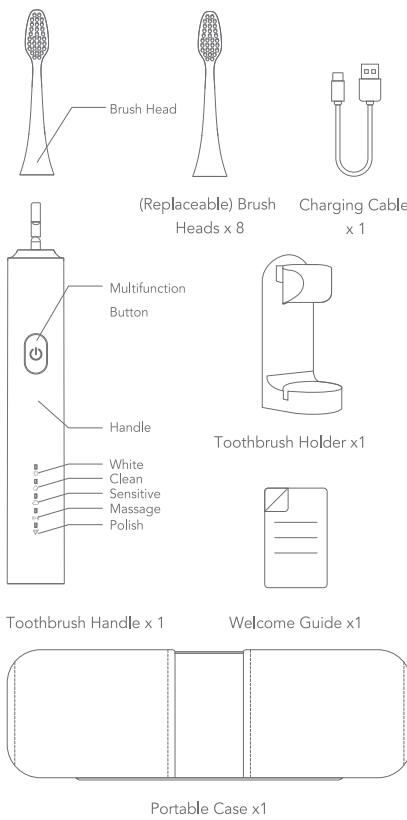
BASIC PRODUCT SPECIFICATIONS

Cover Material: PC + ABS
Battery Capacity: 500mAh
Normal Usage Time: ≥ 60 days
Charging Time: ≤ 6 hours
Rated Working Voltage: DC 3.7V
Rated Working Power: 2W
*oraimo Lab test results

INSTALLATION INSTRUCTIONS

USAGE INSTRUCTIONS

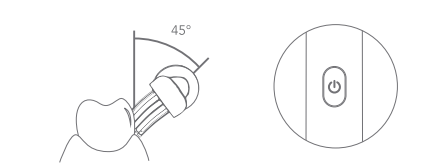Oral Cavity

- Wet the bristles and apply a suitable amount of toothpaste. To prevent splattering, position the toothbrush on your teeth before turning it on.
- Angle the bristles slightly towards the gum line and position them against your teeth, hold the brush at a 45° angle to your teeth and gums.
- Press the power switch and gently glide the brush head across your teeth using small, vertical motions. For each section (as illustrated), brush the front, back, and chewing surfaces for 10 seconds each.
- After brushing, remove the toothbrush from your mouth, rinse it thoroughly, and store it in a well-ventilated, dry location.

MAIN SPECIFICATIONS

Brushing Mode

Gear Mode	Frequency/Speed	Description	Product Status
Gear 1	32000 times/minute	Strong gear	The first white indicator light is on
Gear 2	32000 times/minute	Mid-range	The second white indicator light is on
Gear 3	32000 times/minute	Weak gear	The third white indicator light is on
Gear 4	32000 times/minute	Massage mode	The fourth white indicator light is on
Gear 5	36000-47000 times/minute	Variable speed	The dynamic frequency changes rapidly, and the fifth white indicator light is on

Battery Indicator

Battery Level	Indicator Light After Shutdown	Product Status
Less than 10%	No indicator light	Product not working
10% - 20%	One gear light	Product in low power mode
20% - 40%	Two gear lights	Product in medium power mode
40% - 60%	Three gear lights	Product in medium-high power mode
60% - 80%	Four gear lights	Product in high power mode
Greater than 80%	Five gear lights	Product fully charged

Charge Indicator

Battery Level	Indicator Light Status
Less than 10%	Charging indicator light not on
10% - 90%	Gear lights gradually light up from bottom to top, indicating a cycle
90% - 100%	Charging indicator light not on

TROUBLESHOOTING

Trouble	Solutions
The product does not work.	Check if the battery has power .
The vibrational frequency has decreased.	The battery may be low. Please charge it in time.
The brush head has less cleaning effect.	Replace the toothbrush head at least once every three months.

IMPORTANT SAFETY INSTRUCTIONS

- Do not open, repair, or modify the product without authorization.
- Ensure toddlers or individuals unable to operate the device are supervised by an adult.
- Avoid submerging the product in water for extended periods.
- Remove the battery if the product will not be used for a long time.
- Do not exceed the product's intended functional use.
- This product is designed for household use to clean teeth, gums, and tongue only. Do not use it for other purposes.
- Use the vibrational strength provided by the toothbrush; do not apply excessive pressure on the teeth or gums.
- Keep the product away from high temperatures and direct sunlight.
- Avoid striking or applying excessive force to the product.

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