



WELCOME GUIDE

Sonicare | OPC-ET200N I

*Please read the instructions carefully before using the product and keep the user manual in a safe place for future reference.

SCAN & DOWNLOAD



oraimo store
Easily Smart, Simply Best



Authentic Products



Fast Delivery




Best Price



Prompt Service


PACKAGE CONTENTS & PRODUCT OVERVIEW

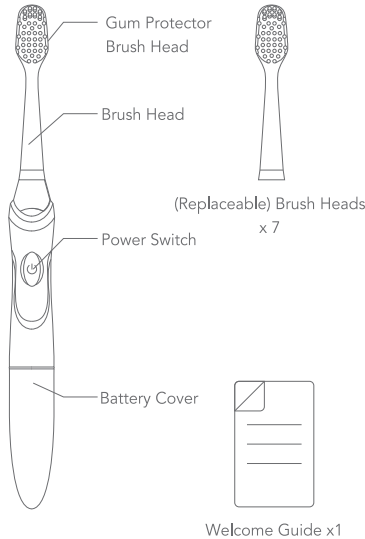


BASIC PRODUCT SPECIFICATIONS

Cover Material: PC+ABS
Rated Voltage: 1.5V
Rated Power: 0.3W
Battery Type: 1 x AAA Battery
Usage Time: ≥ 120 Days
Vibrate Speed: 20,000 Times / Min
Size: 216.5 x 19.8 x 22.2mm
*oraimo Lab test results

INSTALLATION INSTRUCTIONS

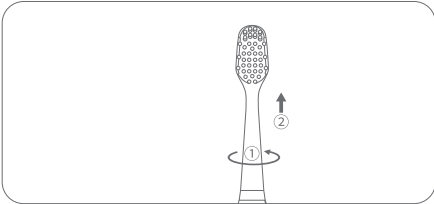

① Step 1: Turn brush head ① and pull up ②.

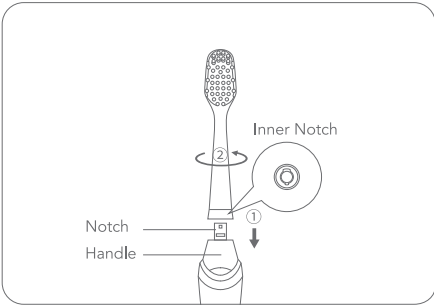

② Step 2: Align the new brush head ① and turn it ②.

Note:

Ensure the notch on the neck aligns with the inner notches of the toothbrush head to lock it in place.

USING INSTRUCTIONS

- Wet the bristles and apply a suitable amount of toothpaste. To prevent splattering, position the toothbrush on your teeth before turning it on.
- Angle the bristles slightly towards the gum line and position them against your teeth.
- Press the power switch and gently glide the brush head across your teeth using small, vertical motions. For each section, brush the front, back, and chewing surfaces for 10 seconds each.
- After brushing, remove the toothbrush from your mouth, rinse it thoroughly, and store it in a well-ventilated, dry location.

TROUBLESHOOTING

| Trouble | Solutions |
|--|--|
| The product does not work. | Ensure an AAA battery is inserted. Check if the battery has power and is correctly installed (positive and negative poles aligned properly). |
| The vibrational frequency has decreased. | The battery may be low. Replace it with a new one. |
| The brush head has less cleaning effect. | Replace the toothbrush head at least once every three months. |

IMPORTANT SAFETY INSTRUCTIONS

- Do not open, repair, or modify the product without authorization.
- Ensure toddlers or individuals unable to operate the device are supervised by an adult.
- Avoid submerging the product in water for extended periods.
- Remove the battery if the product will not be used for a long time.
- Do not exceed the product's intended functional use.
- This product is designed for household use to clean teeth, gums, and tongue only. Do not use it for other purposes.
- Use the vibrational strength provided by the toothbrush; do not apply excessive pressure on the teeth or gums.
- Keep the product away from high temperatures and direct sunlight.
- Avoid striking or applying excessive force to the product.