



## WELCOME GUIDE


FlossJet Pro | OPC-CWF201N |

\*Please read the instructions carefully before using the product and keep the user manual in a safe place for future reference.

### IMPORTANT SAFETY INSTRUCTION

1. Do not disassemble the product or change the parts of the product to avoid damage on the product.
2. Do not use above 50°C hot water into the water tank for cleaning mouth to avoid scald.
3. Do not put the product into the fire or heating and do not charge, use or leave it at the is high temperature.
4. Do not use this product for any other purpose other than oral cleaning.
5. Do not push the nozzle hardly on the gums or teeth to avoid hurts.
6. Do not immerse this product into water for a one time to avoid failure or cause fire, electric shock and other dangers.
7. Do not damage, destroy or modify the charging cable and do not stretch, bend or twist it. Do not press heavy objects on the charging cable or get it caught dipped between objects, and keep the the product and cable from hot objects. Do not store the product into in a hot or humid place and avoid direct sunlight.
8. Do not charge at the product a hot or humid place, such as the bathroom or toilet.
9. If not in use of the product for a long time, please unplug the charging cable to avoid fire or danger.
10. Do not use the product without a nozzle.
11. This product is not suitable for children and infants unless they are supervised or guided by their safety responsible supervisor in the use of this product.
12. Keep out of reach for children or infants. Do not use it as a toy.
13. People who have periodontal disease or have had oral surgery in the past 2 months, please consult your dentist before using.
14. Do not add mouthwash into the water tank which will cause corrosion and the tank will be adhered to the device body can result in tank opening failed.


### PACKAGE CONTENTS & PRODUCT OVERVIEW

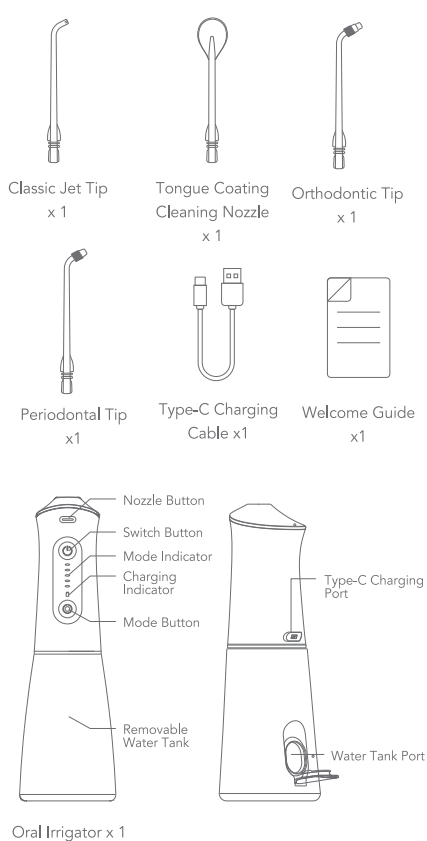

### BASIC PRODUCT SPECIFICATIONS

Material: ABS/PC/PP  
Product Size: 83.53\*65.86\*222.04mm  
Rated Power: 8W  
DC Input: 5V 0.6A (Type-C)  
Charging Time: 6 Hours  
Battery: 2000mAh  
\*oraimo Lab test results

### USAGE INSTRUCTIONS

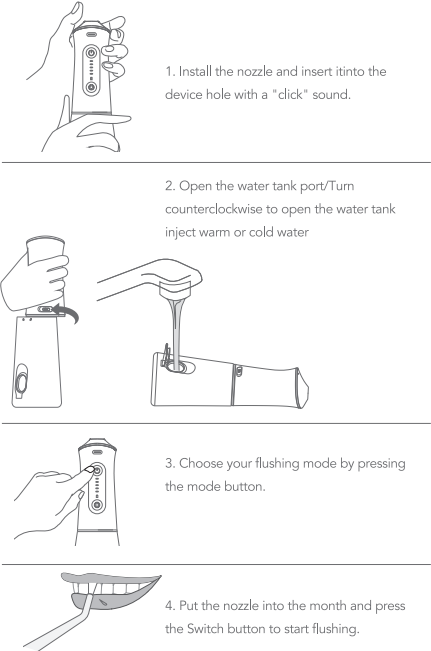

5. Open the tank cover for regular cleaning and keep the tank clean.

### MODE INTRODUCTION

| Mode Indicator<br>(from bottom to top) | Indicator Status                                        |
|----------------------------------------|---------------------------------------------------------|
| Pulse Mode                             | The first indicator light turns on.                     |
| Sensitive Mode                         | The second indicator light turns on.                    |
| Medium Mode                            | The second and third indicator lights turn on.          |
| Normal Mode                            | The second, third, and fourth indicator lights turn on. |
| Intense Mode                           | The second to fifth indicator lights turn on.           |

### CHARGING THE MACHINE

| Status Indicator          | Battery Level |
|---------------------------|---------------|
| ● Green light on (steady) | Fully charged |
| ⚡ Red light flashing      | Being charged |
| ● Red light on (steady)   | Low           |

### CLEANING

#### Cleaning

It is recommended to wash thoroughly after each use and keep it dry before storage to avoid the breeding of bacteria.



1. Close the Oral irrigator.



2. You can use a mixture of vinegar and water to clean the inside of the nozzle and tank, and then wipe off the water stains with a dry cloth.



3. Install the water tank and turn it to the lock mark.

### Nozzle

1. It is recommended to replace the nozzle every three months for optimal performance.
2. If the nozzle becomes deformed, please replace it promptly to ensure proper functionality.

### Suction hose

1. Do not bend, twist or pull the suction hose.

### Body

1. Ensure the charging port is securely closed before cleaning the product body to prevent water ingress.
2. For routine cleaning, use a soft cloth moistened with mild detergent to wipe the surface, then dry with a clean cloth.
3. If the unit is heavily soiled, it can be rinsed with water. Do not use solvents, gasoline, alcohol, or other flammable or corrosive liquids to clean the product.
4. Do not use hot water above 50°C when cleaning.

### Water tank cap

1. Regularly inspect the water tank cap to ensure it is clean and free from residue.
2. Failure to do so may lead to water leakage, bacterial growth, and potential hygiene risks.

### Water Tank

1. Always fill the tank with clean water before each use. After use, it is recommended to remove the tank and clean the inner walls thoroughly to minimize bacterial growth and ensure hygiene.
2. If the product will not be used for an extended period, ensure the water tank is emptied and completely dried before storage